



HAMILTON-G5

The modular high-end ventilation solution

HAMILTON
MEDICAL



We live for ventilation technology

We live for ventilation technology. Technology that helps caregivers improve the lives of their critically ill patients. We believe that innovation is essential to meet the demands of critical care. To us, innovation is about realizing visionary new ideas and continuously improving existing products, always maintaining the focus on safe, individualized ventilation, as well as ease of use.

We learn from our customers and from multi-disciplinary experts. And we invest in long-term research and development. We develop Intelligent Ventilation solutions: devices and consumables for the ventilation of all critically ill patients – from neonates to adults.

A handwritten signature in blue ink that reads "Jens Hallek". The signature is fluid and cursive.

Jens Hallek
CEO
Hamilton Medical AG

A handwritten signature in blue ink that reads "Bob Hamilton". The signature is fluid and cursive.

Bob Hamilton
CEO
Hamilton Medical, Inc.

Meet the HAMILTON-G5

The HAMILTON-G5 is Hamilton Medical's most modular high-end mechanical ventilator. A wide range of standard features and options allows you to tailor the HAMILTON-G5 to your needs.

- ✓ Real-time patient synchronization with IntelliSync+
- ✓ P/V Tool® Pro for lung assessment and recruitment
- ✓ Transpulmonary pressure measurement
- ✓ High flow oxygen therapy
- ✓ Integrated IntelliCuff pressure controller
- ✓ Remote access to humidifier controls and status
- ✓ Adult, pediatric, and neonatal ventilation
- ✓ Adjustable O2 enrichment
- ✓ Adjustable volume limitation for neonatal patients



Flexible device configuration

Customized solution

You can configure the HAMILTON-G5 in several different ways to customize the device for your environment. If using a trolley, you can mount the monitor either on top or in front of the ventilation unit. If the unit is positioned on a shelf, you can mount the monitor on top or to the side using the side-mount option. You can also adjust the orientation and angle of the monitor by turning and tilting it as required.

Optimal visibility

The 360°-visible alarm lamp on top of the HAMILTON-G5 monitor allows you to identify alarms easily. The optional nurse-call capability provides additional support for optimal alarm detection. The 15-inch touch screen was designed for smooth and fast operation. It also gives you an overview of the patient's current ventilation status at a quick glance and provides a reliable basis for therapy decisions.





Ease of use

In close cooperation with users and ventilation experts, our engineers have designed a user interface that is particularly intuitive. Switching between the HAMILTON-G5 and all other Hamilton Medical ventilators is easy because they are all operated according to the same principles.

The Ventilation Cockpit on the HAMILTON-G5 consolidates the monitoring data and displays it as intuitive graphics. These provide a quick overview of the patient's current ventilation status and provide a reliable basis for therapy decisions.



What I like the most about the HAMILTON-G5s is probably the monitoring parameters and the ability to trend those in real time and up to 72 hours. I have been able to use that on a disease-specific basis and trend data that I couldn't do before.

Craig Jolly, Adult Clinical Education Coordinator
University Medical Center, Lubbock (TX), USA



The Ventilation Cockpit

1 Main monitoring parameters

All of the main monitoring parameters and alarm limits at a glance. The large characters allow you to see them even from a distance.

2 Dynamic Lung

One quick look shows you tidal volume, lung compliance, resistance, and patient efforts/ triggers in real time. The lungs expand and contract in synchrony with the actual breaths.

3 Vent Status

The Vent Status panel displays six parameters related to the patient's dependence on the ventilator. When all values are in the weaning zone, the panel is framed in green, indicating that spontaneous breathing trials or extubation can be considered.

4 Direct access to main controls

Access and adjust the most important controls for the current mode directly on the main display.



Individualized, lung-protective ventilation

The features available on the HAMILTON-G5 help you to individualize your patient's ventilation and to implement a lung-protective ventilation strategy.

Adaptive, lung-protective ventilation with ASV

- ✓ Supports the earliest possible spontaneous breathing by the patient^{1, 2}
- ✓ Shortens the ventilation time in various patient groups^{1, 2}

Lung assessment and recruitment with the PV Tool Pro

- ✓ Hysteresis of the pressure/volume curve can be used for assessing the recruitability of the lung at the bedside³

Synchronization based on waveform analysis with IntelliSync+

- ✓ Waveform analysis is a reliable, accurate, and reproducible method for assessing patient-ventilator interaction⁴
- ✓ In terms of cycling, IntelliSync+ performs at least as well as ETS optimized by clinicians⁵

Automatic cuff pressure control with IntelliCuff

- ✓ Continuous cuff pressure control can decrease microaspiration and VAP^{6, 7}

Transpulmonary pressure measurement

- ✓ PEEP set based on transpulmonary pressure can improve compliance and oxygenation in ARDS patients⁸
- ✓ Transpulmonary pressure measurement can avoid the use of ECMO in the most severe patients⁹

¹ Kirakli C. Eur Respir J. 2011 Oct;38(4):774-80

² Chen CW. Respir Care. 2011 Jul;56(7):976-83

³ Demory D. Intensive Care Med. 2008 Nov;34(11):2019-25

⁴ Mojoli F. Intensive Care Med Exp 2016, 4(Suppl 1):A1168

⁵ Mojoli F. Intensive Care Med Exp 2016, 4(Suppl 1):A1164

⁶ Lorente L. Critical Care. 2014;18(2):R77

⁷ Nseir S. American Journal of Respiratory and Critical Care Medicine. 2011;184(9):1041-1047

⁸ Talmor D. N Engl J Med. 2008 Nov 13;359(20):2095-104

⁹ Grasso S. Intensive Care Med. 2012 Mar;38(3):395-403



Adaptive Support Ventilation (ASV)

adjusts respiratory rate, tidal volume, and inspiratory pressure continuously, depending on the patient's lung mechanics and effort. ASV adapts ventilation breath-by-breath, 24 hours a day, and from intubation to extubation.



P/V Tool Pro for lung assessment and recruitment

helps you assess recruitability and set PEEP based on respiratory mechanics. It also provides a repeatable method for quickly performing recruitment maneuvers.



IntelliSync+ keeps an eye on patient-ventilator synchrony

by continuously analyzing waveform shapes hundreds of times per second. This allows IntelliSync+ to detect patient efforts and cycling immediately, and initiate inspiration and expiration in real time. IntelliSync+ applies to invasive and noninvasive ventilation, regardless of the ventilation mode.



IntelliCuff pressure controller

continuously measures and automatically maintains the user-set cuff pressure of an endotracheal or tracheostomy tube in real time.



Transpulmonary pressure measurement

allows optimization of PEEP, tidal volume, and inspiratory pressure. Use it in combination with the P/V Tool Pro to assess lung recruitability more precisely and perform recruitment maneuvers.

Features and options



State-of-the-art ventilation modes



Adult, pediatric, and neonatal ventilation



Remote access to HAMILTON-H900 controls and status



Integrated control for IntelliCuff pressure controller



Integrated pneumatic and optional Aerogen nebulizer




Heliox therapy



Pulse oximetry (SpO₂ and pulse measurement)



Continuous monitoring of driving pressure

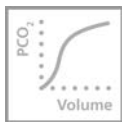

Mainstream (volumetric) and sidestream capnography



Configurable loops and trends



Serial interface for connection to PDMS or patient monitors



Hot-swappable battery backup



High-performance noninvasive ventilation (NIV)



Integrated high flow oxygen therapy

From the ventilation specialist

E-learning

Hamilton Medical College provides free and open e-learning on mechanical ventilation and ventilators.

Join at: www.hamilton-medical.com/elearning.

Universal ventilator consumables

Our accessories and consumables are specially developed for the highest possible patient safety and ease of use. Choose between reusable and disposable parts according to your institutional policies.

Peripheral devices

Our ventilation portfolio includes an active humidifier, the HAMILTON-H900, as well as the automatic cuff pressure controller, IntelliCuff. Both devices may be used with all kinds of mechanical ventilators.





More information and free simulation software:
www.hamilton-G5.com



HAMILTON **MEDICAL**

Hamilton Medical, Inc.

4655 Aircenter Circle, Reno, NV 89502

☎+1 800 426 6331

info@hamiltonmedical.com

www.hamilton-medical.com

ELO201604045.03

Specifications are subject to change without notice. Some features are options. Not all features are available in all markets. INTELLiVENT-ASV is not available in the US. For all proprietary trademarks and third-party trademarks used by Hamilton Medical AG see www.hamilton-medical.com/trademarks. ©2021 Hamilton Medical AG. All rights reserved.

HAMILTON-G5