

The basics of high flow nasal cannula therapy



Synonyms and abbreviations

There is a range of terminology used for high flow nasal cannula (HFNC) therapy: HFNC supportive therapy, high flow therapy (HFT), nasal high flow (NHF), high flow (HF), high flow nasal cannula oxygen therapy, heated humidified high flow therapy (HHHF), high flow oxygen (HFO) therapy, and high flow oxygen therapy (HFOT). Hamilton Medical uses the term **high flow nasal cannula (HFNC) therapy**.

ARF	acute respiratory failure
CO ₂	carbon dioxide
COT	conventional oxygen therapy
EELV	end-expiratory lung volume
ETI	endotracheal intubation
FiO ₂	fraction of inspired oxygen
NIV	noninvasive ventilation
O ₂	oxygen

PaO ₂	partial pressure of oxygen
PEEP	positive end-expiratory pressure
PIF	peak inspiratory flow
RR	respiratory rate
SpO ₂	oxygen saturation
VE	minute volume
WOB	work of breathing



Note: The term HFNC refers to the therapy itself irrespective of the interface used (nasal cannula or tracheostomy connector).

Content overview

Your ventilation expert	4
Introduction	5
What are the benefits?	8
How does it work?	10
What is the clinical relevance?	18
How to choose the right interface?	21
How to adjust the settings?	23
How to monitor treatment?	27
Test yourself	32
Appendix	44
References	50

Your ventilation expert



Kathrin Seeger

ICU nurse

Hamilton Medical

Clinical Application Specialist

Introduction

High flow nasal cannula therapy is a type of respiratory support that delivers heated and humidified gas with a controlled concentration of oxygen to your patients. This therapy has become increasingly popular for hypoxemic respiratory failure.

HFNC is used during early noninvasive management of acute respiratory failure, and has proven to be safe and effective as a noninvasive ventilation method¹.


HFNC has been shown to be potentially useful and efficient in other applications such as major post-operative care, immunocompromised patients, for preoxygenation, or during bronchoscopy².



This e-book focuses on the basics of high flow nasal cannula therapy for your daily clinical practice.

Introduction

The system used to deliver HFNC requires the following components: a gas blender and flow meter, an active humidifier, a heated inspiratory circuit, and a nasal cannula or a tracheostomy connector as an interface. When we talk about high flow nasal cannula, we need to keep in mind the four elements that describe this therapy: flow, oxygen, heat, and humidity.

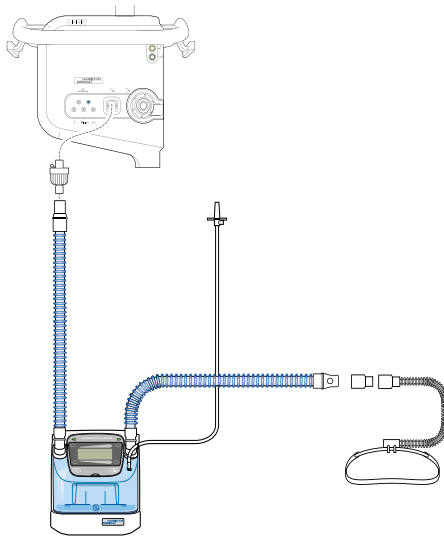

The four elements of HFNC

Flow



Oxygen



Heat



Humidity



Introduction

If we take a closer look at the meaning of each element, the advantages of using HFNC are clear:



Flow

The higher flow rates allow you to set the inspiratory flow to adequately meet or even exceed the inspiratory demand of the patient.



Oxygen

You can set the FiO₂ more accurately within a range of 21% to 100%. However, consider escalating treatment if FiO₂ of greater than 60% is needed to achieve the targeted SpO₂.



Heat

The gas should be heated to the core body temperature of 37°C (except during NIV where the default temperature is 31°C).



Humidity

The breathing gas is saturated close to 100% relative humidity to meet the physiologic demands.



Did you know that the inspiratory demands of adult patients with acute respiratory failure can range from 30 l/min to more than 120 l/min?

What are the benefits?

What are the benefits

We can divide the benefits of HFNC into clinical and physiologic as follows:

Clinical benefits¹

- ✓ Improved patient comfort and compliance with the treatment
- ✓ Better patient tolerance due to the warmth and humidity of inspired air
- ✓ Ease of use

Physiologic benefits¹


- ✓ Improved oxygenation because higher inspiratory flows mean reduced oxygen dilution
- ✓ Reduction of dead space and CO₂ rebreathing compared with low flow oxygen therapy and NIV
- ✓ Reduced **work of breathing** ↻ associated with an improvement in respiratory rate, as well as a reduction in heart rate with a significant improvement in SpO₂ measurements

How does it work?



From low flow to high flow

Before we dive into the working principles of high flow nasal cannula, we need to address the limitations that are typical of low flow oxygen delivery systems. This will help you better understand the beneficial effects of HFNC for your patients.

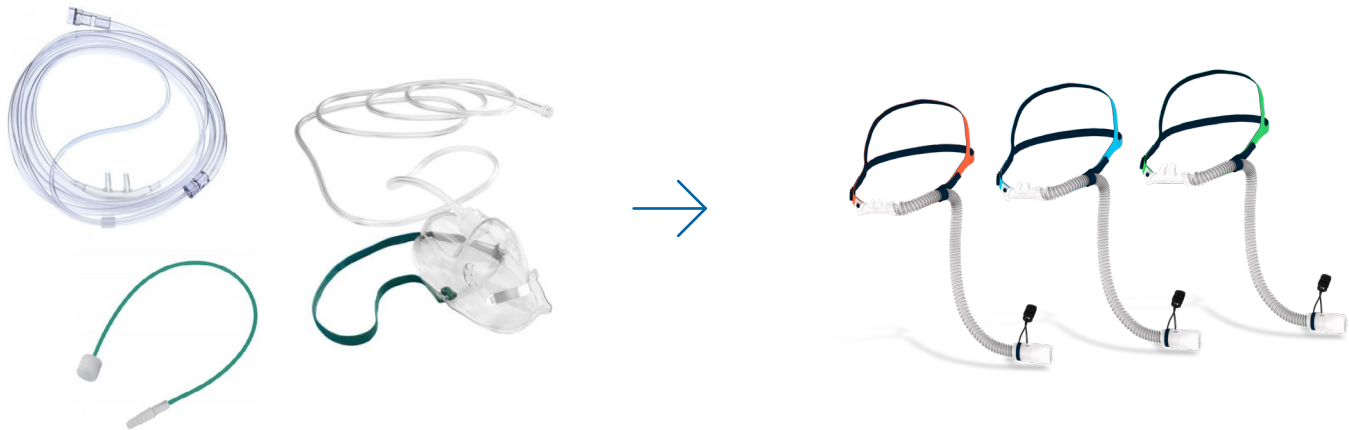

Low flow oxygen therapy

With conventional oxygen delivery devices, the peak inspiratory flow rate in acute respiratory failure patients exceeds the delivered flow of oxygen. Typical low flow systems provide supplemental O₂ directly to the airway at a flow of 8 l/min or less. This O₂ provided by a low flow device is always diluted with ambient air (with 21% FiO₂); the result is a low and variable FiO₂.

This variable FiO₂ may be associated with one or more of the following: the interface used to deliver oxygen, the respiratory rate, the peak inspiratory flow and how the patient breathes (through the mouth or nose).



Low flow nasal cannulas

Only 22% to 40% of oxygen at flow rates of up to 6 l/min in adults



Low flow oxygen masks

Only 35% to 50% of oxygen at flow rates from 5 to 10 l/min



The use of humidification is recommended when oxygen is supplied through a nasal cannula at flow rates of more than 4 l/min¹.


High flow nasal cannula

- In contrast to low flow oxygen therapy, the gas flow rate and FiO_2 in HFNC can be adjusted independently of one another, depending on the patient's inspiratory demands
- With high flow nasal cannula systems, the FiO_2 is delivered more accurately and can be set from 21% to 100%
- As mentioned already, the higher flow rates are able to meet or even exceed the patient's peak inspiratory flow rate
- The blend of air and oxygen delivered to the patient is fully heated and humidified by the humidification chamber and the single limb heated breathing circuit. This breathing circuit contains heater wires within the tubing wall that minimize condensation
- The therapy is delivered directly into the patient's nostrils through a nasal cannula



Low flow vs. high flow nasal cannula therapy

During normal inspiration, the PIF demand is 30 to 40 l/min. With high flow nasal cannula therapy, the administered FiO_2 would equal the inhaled FiO_2 , while in the case of a low flow system, the inhaled FiO_2 would be less than the administered FiO_2 because it is diluted with ambient air³.

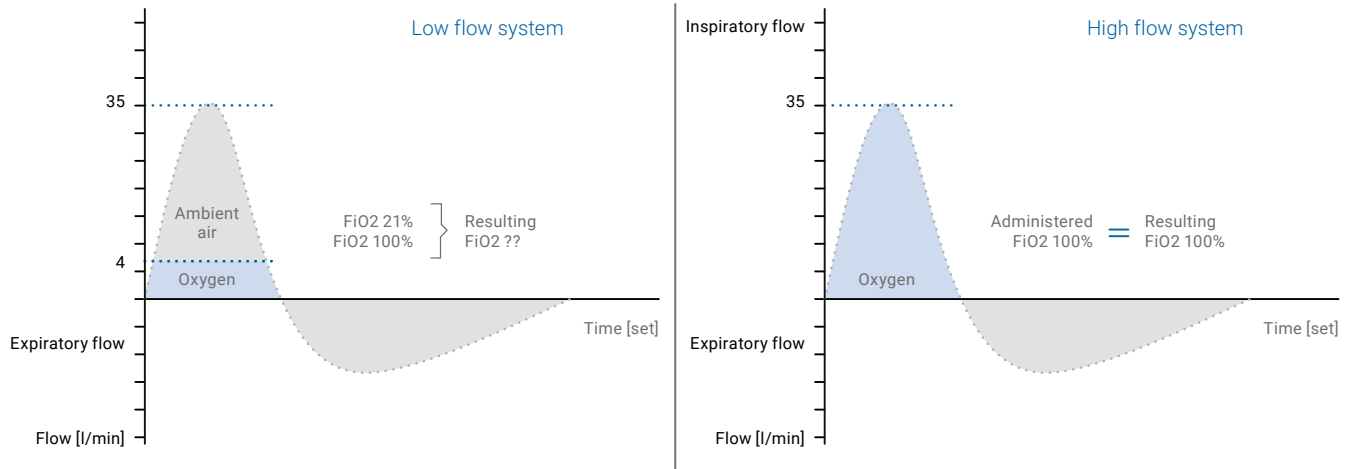


Low flow vs. high flow nasal cannula therapy

For higher flows, adequate humidification is fundamental because it improves the tolerance and compliance of the treatment. HFNC systems require an active heated humidifier to achieve optimal body temperature (37°C in adults and 34°C to 37°C in pediatric patients) with a humidity output of higher than 33 mg H₂O/l and close to 100% of relative humidity.


These conditions provided by HFNC promote:

- ✓ Improved lung compliance¹
- ✓ Improved gas distribution in the lungs¹
- ✓ Improved mucociliary clearance¹
- ✓ Decreased airway resistance¹
- ✓ Increased patient comfort¹

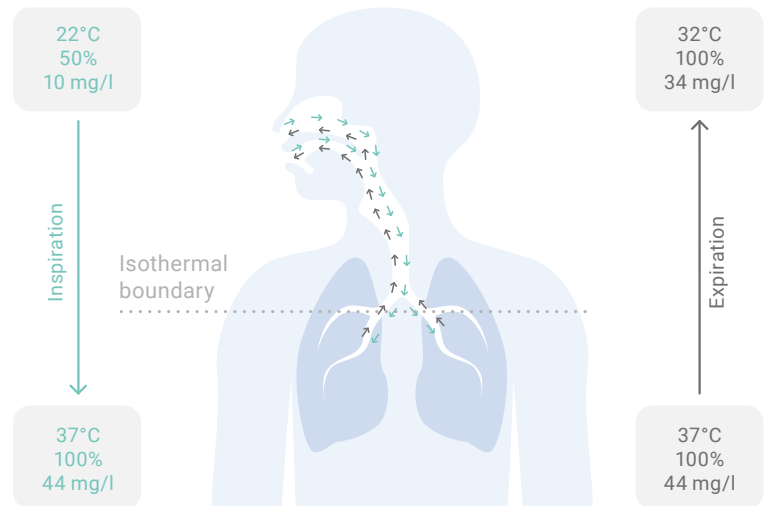

When to consider high flow nasal therapy

HFNC is considered the first-line therapy for patients with acute hypoxemic respiratory failure and can be part of a treatment strategy, whether during escalation or weaning.

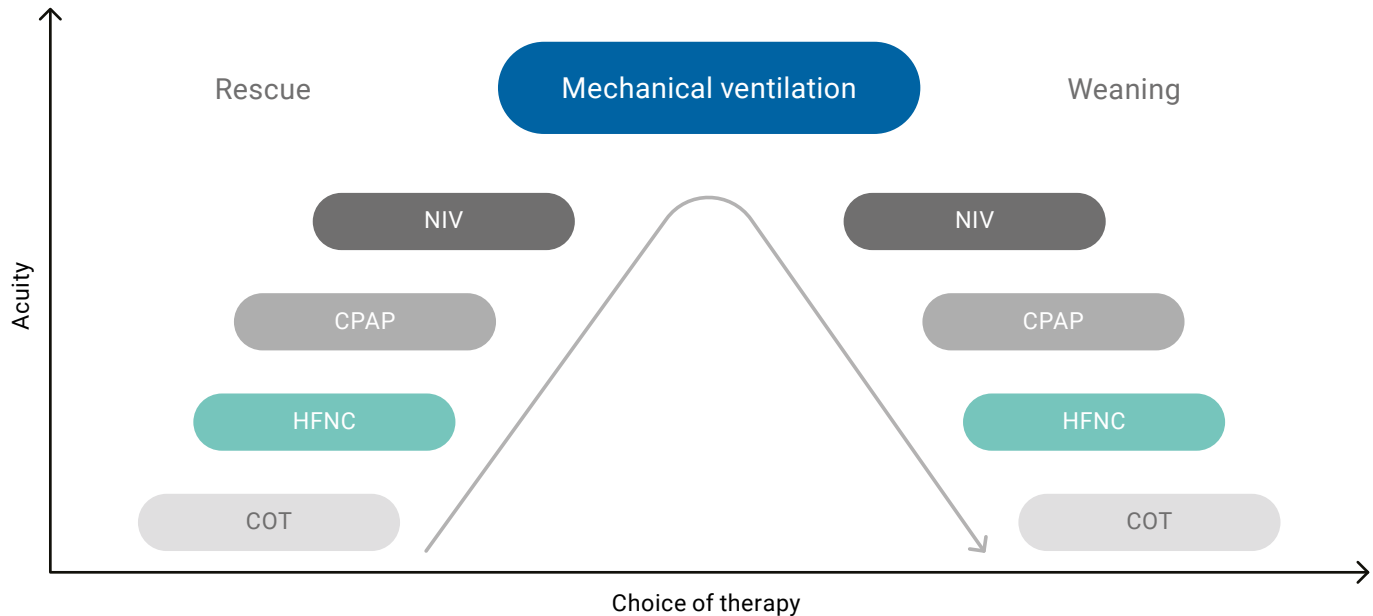

Recommendations for HFNC^{4,5}

Hypoxemic
respiratory failure



Strong
recommendation

- ✓ Lower intubation rates
- ✓ Reduction in escalation of respiratory support
- ✓ Cost savings in terms of equipment and intubations avoided

Post-extubation



Conditional
recommendation

- ✓ Lower intubation rates (compared with COT)
- ✓ Reduction in **post-extubation respiratory** ↻ failure

Post-operative HFNC in
high risk and/or obese
patients following cardiac
or thoracic surgery



Conditional
recommendation

- ✓ Lower reintubation rate (compared with COT)
- ✓ Reduction in escalation of respiratory support (compared with COT)

Peri-intubation period



No recommendation

- ✓ Patients who are already receiving high flow should continue with HFNC during intubation

What is the clinical relevance?



Physiologic effects of high flow

Negative end-expiratory transpulmonary pressure can cause ventilator-induced lung injury due to atelectrauma. Studies⁶ have shown that HFNC has multiple physiologic effects ranging from improved lung volumes to decreased work of breathing. The main ones are as follows:


- ✓ High flow generates a washout effect in the upper airways, promoting anatomic dead-space clearance and CO₂ removal
- ✓ Increases the end-expiratory lung volume and PaO₂/FiO₂ ratio
- ✓ Optimal strategy for administering oxygen to hypoxic critically ill patients with high respiratory demand
- ✓ Provision of low-level PEEP
- ✓ The therapy is delivered directly into the patient's nostrils through the nasal cannula



The physiologic effects of HFNC can be summarized as improved oxygenation, reduced work of breathing, improved lung protection and better comfort for the patient.

The basics of high flow nasal cannula therapy

20

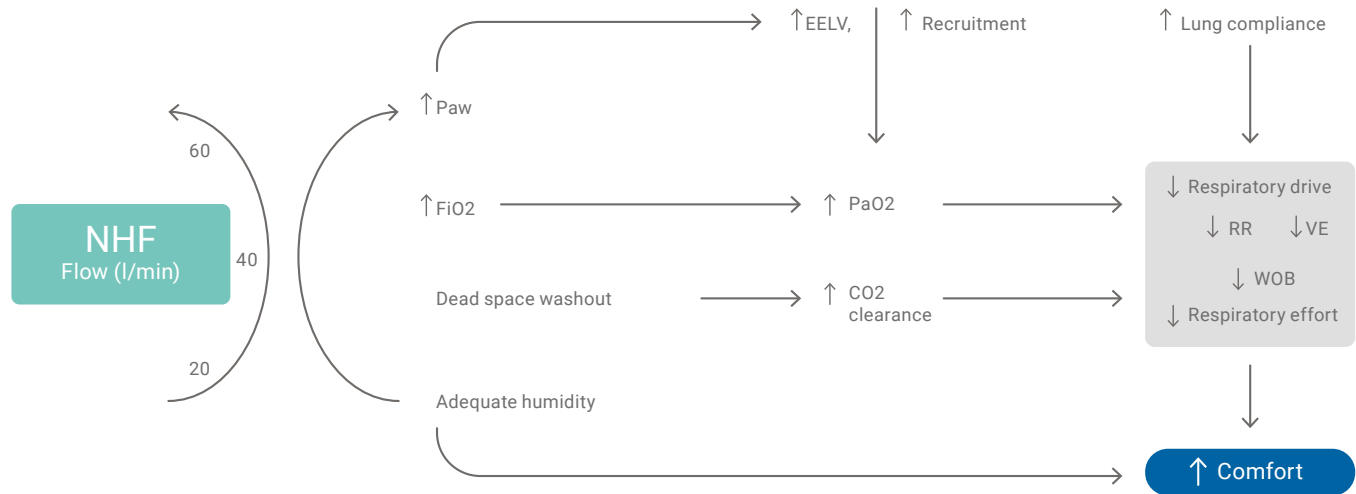


How to choose the right interface?



Choosing the right cannula size

The nasal cannula interface is an essential component for delivery of high flow nasal cannula therapy. In order to maintain an effective flush of CO₂, it is important that the cannula does not occlude more than ~50% of the nares.

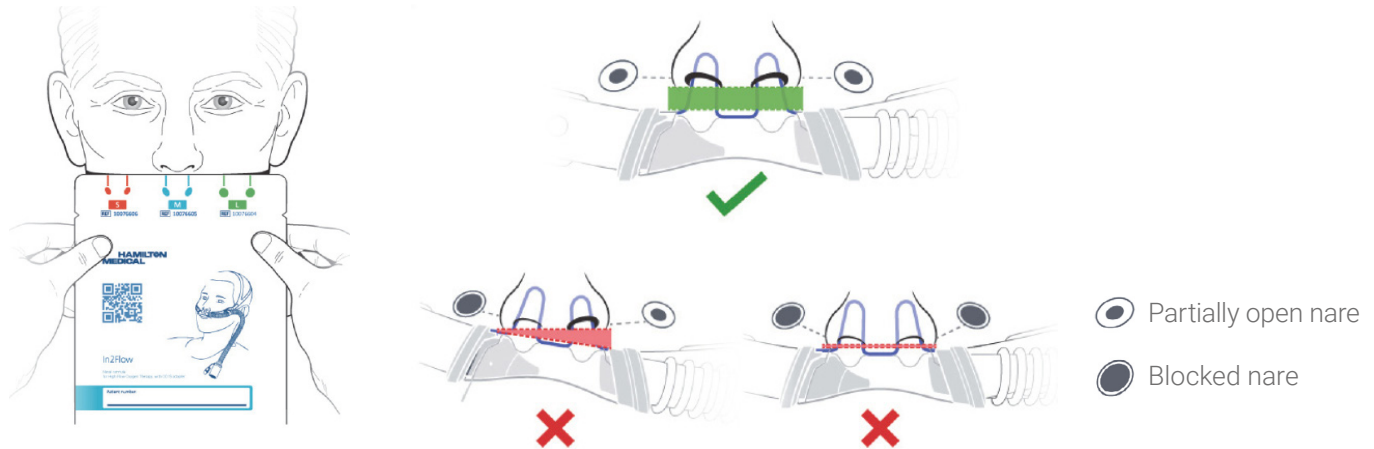

Note: Airway pressure increases progressively with both increasing flow rate and nasal prong-to-nares ratio.

How to adjust the settings?



Setting the flow in adult patients⁷

- ✓ Set the flow [↗](#) initially to between 20 and 35 l/min
- ✓ The flow can be increased progressively in steps of 5 to 10 l/min if the respiratory rate fails to improve
- ✓ Increasing flows from 15 to 45 l/min triples the reduction in anatomic dead space
- ✓ To prevent intubation, choose the highest flow tolerated by the patient

These are only general recommendations, and every patient should be treated based on their specific medical condition.



Setting the flow in pediatric patients⁸

Flow rates that exceed inspiratory demand can be set in patients under 24 months who tolerate flows of 1 to 2 l/kg/min (up to 20 l/min).

Age	Body weight (kg)	Flow range (l/min)
≤ 1 month	< 4	5–8
1–12 months	4–10	8–20
1–6 years	10–20	12–25
6–12 years	20–40	20–30
12–18 years	> 40	25–50

These are only general recommendations, and every patient should be treated based on their specific medical condition.

Setting the oxygen in all patient groups

The FiO₂ is the concentration of oxygen in the gas mixture delivered to the patient. In HFNC, the FiO₂ should be set as follows:

- ✓ You can set the FiO₂ from 21% up to 100%
- ✓ Titrate the FiO₂ to achieve the desired SpO₂ (target ranges of 92%–96% for most patients and 88%–92% for patients with chronic respiratory disease)⁷

These are only general recommendations, and every patient should be treated based on their specific medical condition.

How to monitor treatment?



Monitoring the effectiveness of HFNC

The ROX index is defined as the ratio of oxygen saturation as measured by SpO₂/FiO₂ to the respiratory rate. In [patients with hypoxemic respiratory failure](#) ⁹, the ROX index can be used after HFNC therapy has been initiated to help identify those who are at high risk of intubation⁹.

ROX index ≥ 4.88 after two hours of treatment, indicates a high probability that intubation will not be necessary.

ROX index < 3.85 indicates a higher risk of treatment failure.

ROX index

$$\text{ROX} = (\text{SpO}_2/\text{FiO}_2)/\text{respiratory rate}$$

ROX ≥ 4.88 success

SpO ₂ :	94%
FiO ₂ :	0.6
SpO ₂ /FiO ₂ :	157
Rate :	25

ROX: 6.27

ROX < 3.85 consider intubation

SpO ₂ :	92%
FiO ₂ :	0.8
SpO ₂ /FiO ₂ :	115
Rate :	35

ROX: 3.29

Recommendations for measuring the ROX index⁹




SpO2 monitoring

SpO2 is one of the most important parameters to monitor during HFNC. It indicates whether FiO2 needs to be re-adjusted and is useful during FiO2 titration.

How to use [SpO2* measurement](#) on your Hamilton Medical ventilator:

- ✓ Connect the sensor to the device and the patient
- ✓ Select the System tab
- ✓ Activate the SpO2 sensor

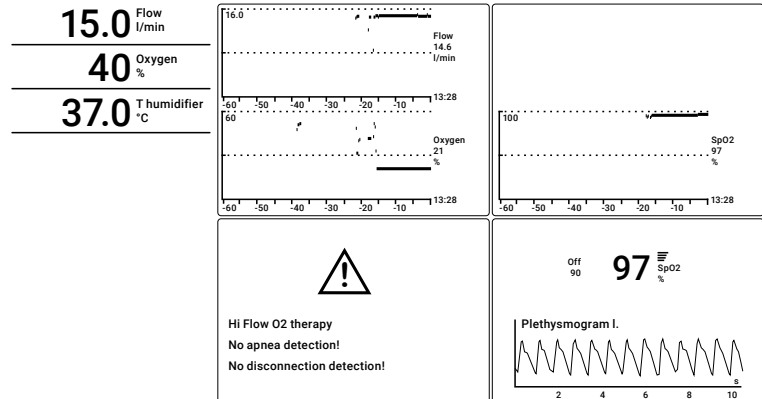

*If the SpO2 option is not installed, use external patient monitoring

SpO2 monitoring

The monitoring shows the values of flow and oxygen over time so you can evaluate the therapy's progress.


How to monitor the progress of the therapy with your Hamilton Medical ventilator:

- ✓ Main monitoring parameters
- ✓ Control flow
- ✓ Oxygen

Trending graphs:

- ✓ Select either Control Flow or Oxygen

35 Control Flow <small>l/min</small>	Flow: 15 l/min Oxygen: 22%			
22 Oxygen <small>%</small>	26 18			
97 SpO2 <small>%</small>	Info	Tests&Calib	Sensors	Settings
	On/Off SpO2	<input checked="" type="radio"/> O2 sensor <input type="radio"/> CO2 sensor <input checked="" type="radio"/> SpO2 sensor active Masimo SpO2 cable: OK SpO2 sensor: OK		

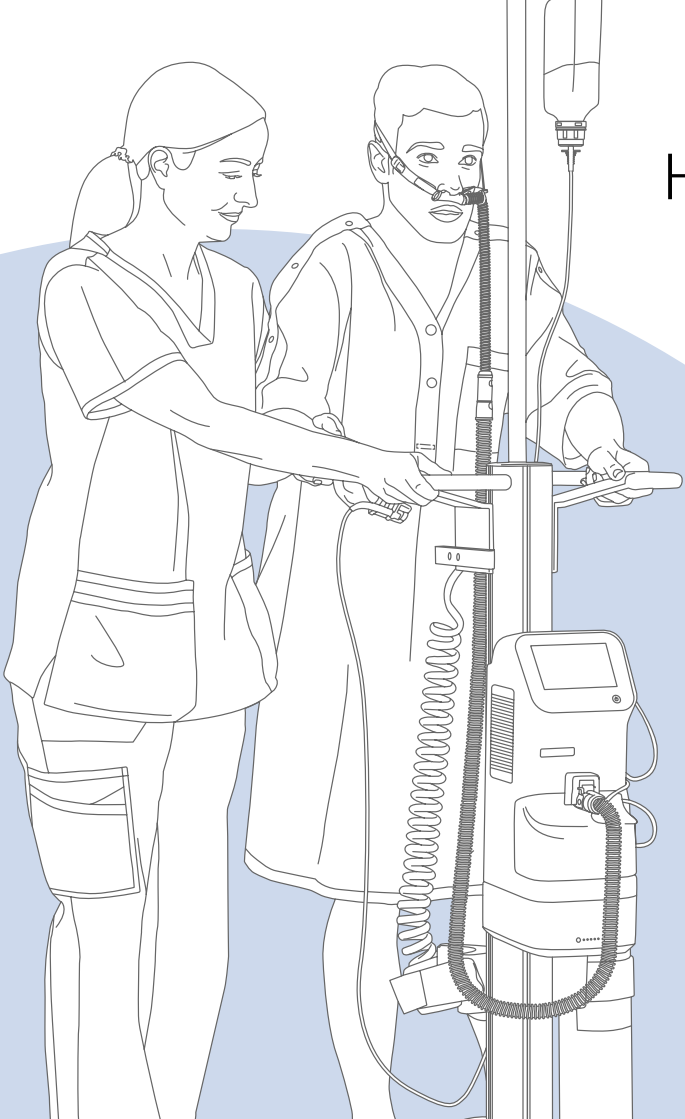

HAMILTON-HF90. Advanced high flow oxygen therapy

www.hamilton-medical.com/HAMILTON-HF90

Learn more

Test yourself



Test yourself

Now it is time to put your knowledge to the test. On the following pages you can find several questions about high flow nasal cannula therapy and the concepts mentioned in this e-book.

For each question there are three possible answers, only one of which is correct.


You can check your answers on page 37 and page 43 (for the clinical case).

Good luck!

A ☐

B ☒


C ☐

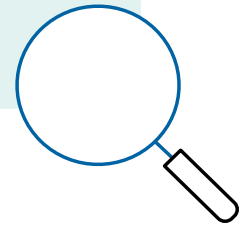

Question 1

When should we consider using high flow nasal cannula therapy?

- a) Preoperative patients
- b) Patients with acute hypoxemic respiratory failure
- c) Patients with chronic hypercapnia

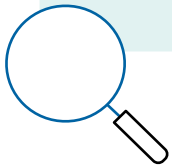

Question 2


Which of the following is a benefit of high flow nasal cannula therapy?

- a) Increased work of breathing
- b) Decreased end-expiratory lung volume and $\text{PaO}_2/\text{FiO}_2$ ratio
- c) Washout effect in the upper airways promoting anatomic dead-space clearance and CO_2 removal

Question 3

What should you set the flow rate to when starting high flow nasal cannula therapy?

- a) 20 to 35 l/min
- b) 15 to 20 l/min
- c) 60 l/min

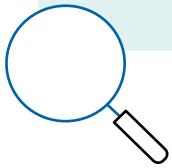

Question 4


How is the ROX index calculated?

- a) $(RR/SpO_2)/FiO_2$
- b) $(SpO_2/FiO_2)/RR$
- c) $(FiO_2/SpO_2)/HF$

Question 5

Which ROX index values indicate an increased risk of treatment failure?

- a) ROX index < 4.88
- b) ROX index < 5.35
- c) ROX index < 3.85

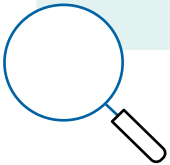

Solutions

Question 1

b) Patients with acute hypoxemic respiratory failure

Question 2

c) Washout effect in the upper airways promoting anatomic dead-space clearance and CO₂ removal

Question 3

a) 20 to 35 l/min

Question 4

b) (SpO₂/FiO₂)/RR

Question 5

c) ROX index < 3.85

Clinical case

Scenario: Emergency transfer of a patient from the normal ward to the ICU due to acute respiratory failure. Following a partial liver resection three days ago, the patient now has rising inflammatory marker levels and a suspected septic event. The patient arrives on the ward with tachypnea and an SpO₂ of 84%. He is hemodynamically stable with a heart rate of 118 bpm and blood pressure of 90/60 mmHg. He is administered 10 l/min of oxygen through an oxygen mask.

Current blood gases:

PaO₂: 72 mmHg (9.6 kPa)


PaCO₂: 52 mmHg (6.93 kPa)

pH: 7.34

Clinical case – Question 1

What would be your first-line therapy to improve the patient's respiratory situation?


- a) NIV
- b) Intubation
- c) HFNC

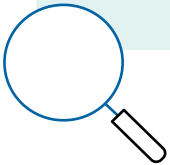

Clinical case – Question 2

For HFNC therapy in this patient, how high would you set the flow?

- a) Between 25 and 40 l/min
- b) Titrate the flow to as high as the patient can tolerate
- c) Flow rate above 60 l/min would be appropriate

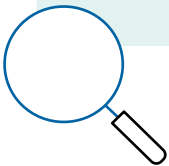

Clinical case – Question 3

At what oxygen level would you start the therapy?


- a) 21%
- b) 50% to 80%; after reaching the desired SpO₂, oxygen is lowered to maintain the SpO₂
- c) Always start with 100%



Clinical case – Question 4

After 30 min of HFNC therapy, RR decreases and SpO₂ is stable at 90%. You decide to continue with the therapy. After 2 hours, you measure the ROX index. The result is 5.93 (SpO₂ - 89% , FiO₂ - 0.6, RR - 25 bpm). How would you proceed now?

- a) Continue therapy without changes
- b) Consider ETI
- c) Change to NIV

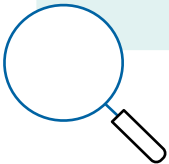

Solutions - Clinical case

Question 1

c) HFNC

Question 2

b) Titrate the flow to as high as the patient can tolerate

Question 3

b) 50% to 80%; after reaching the desired SpO₂, oxygen is lowered to maintain the SpO₂

Question 4

a) Continue therapy without changes

Appendix



HFNC and Hamilton Medical



Hamilton Medical ventilator:

- ✓ With HiFlowO2 mode as an easy-to-add software option
- ✓ Remote access to HAMILTON-H900 controls and status (except HAMILTON-C3)

HAMILTON-H900 humidifier:

- ✓ Active humidification with dedicated high flow mode
- ✓ Preassembled breathing circuit sets with integrated temperature probe

High flow interfaces:

- ✓ Interfaces for adult and pediatric* patients
- ✓ Nuflow nasal cannulas for neonatal patients

*Pediatric patients > 2 years old



HFNC and Hamilton Medical



The [HAMILTON-HF90](#) [↗](#) is an all-in-one solution for high flow oxygen therapy, providing caregivers with the tools and flexibility needed to tailor treatment to each patient's unique needs. Whether in the emergency room, intermediate care, or the ICU, the HAMILTON-HF90 is designed to deliver high flow oxygen therapy whenever required during transfer and mobilization of patients within healthcare facilities.

It combines comprehensive functionality with an intuitive design to support optimal outcomes, enhance patient comfort, and enable patient mobility.



HFNC and Hamilton Medical

Device	 Flow range (l/min)	 Flow range (l/min)
HAMILTON-C1/T1/MR1 ↗	2–30 US market: limited to 15 l/min	2–100 US market: limited to 60 l/min
HAMILTON-C3 ↗	1–12	2–80
HAMILTON-C6 ↗	2–30	2–100
HAMILTON-HF90 ↗	2–30 US market: 2–30	4–80 US market: 4–60

 Neonatal patients
  Adult and pediatric* patients

*Pediatric patients > 2 years old

HAMILTON-H900 humidifier

The [HAMILTON-H900 humidifier](#) was developed with the focus on ease of use and patient safety. This means you can concentrate on other important aspects of patient care.

- ✓ Integrated temperature probe
- ✓ Wall-heated, all-in-one breathing circuit sets
- ✓ Adjustable temperature and humidity settings
- ✓ Remote access to humidifier controls and status*
- ✓ Predefined settings for high flow nasal cannula therapy



*Not available in all markets.

References

1. Carlucci A, Maggiore SM, editors. High Flow Nasal Cannula. Cham: Springer International Publishing; 2021.
2. Ischaki E, Pantazopoulos I, et al. Nasal high flow therapy: a novel treatment rather than a more expensive oxygen device. *European Respiratory Review* 2017 26: 170028.
3. Masclans J, Roca O, et. al. High-flow Oxygen in Acute Respiratory Failure. *Clinical Pulmonary Medicine*. 2012 May;19(3):127.
4. Rochweg B, Einav S, Chaudhuri D, Mancebo J, Mauri T, Helviz Y, et al. The role for high flow nasal cannula as a respiratory support strategy in adults: a clinical practice guideline. *Intensive Care Medicine*. 2020 Dec 17;46(12):2226–37.
5. Oczkowski S, Ergon B, Bos L, Chatwin M, Ferrer M, Gregoretti C, et al. ERS clinical practice guidelines: high-flow nasal cannula in acute respiratory failure. *European Respiratory Journal*. 2022 Apr;59(4):2101574.
6. Ricard JD, Roca O, Lemiale V, Corley A, Braunlich J, Jones P, et al. Use of nasal high flow oxygen during acute respiratory failure. *Intensive Care Medicine*. 2020 Dec 8;46(12):2238–47.
7. Drake MG. High-Flow Nasal Cannula Oxygen in Adults: An Evidence-based Assessment. *Ann Am Thorac Soc*. 2018 Feb;15(2):145–55.
8. Kwon JW. High-flow nasal cannula oxygen therapy in children: a clinical review. *Clinical and Experimental Pediatrics*. 2020 Jan 15;63(1):3–7.
9. Roca O, Messika J, Caralt B, García-de-Acila M, Sztrymf B, Ricard JD, et al. Predicting success of high-flow nasal cannula in pneumonia patients with hypoxemic respiratory failure: The utility of the ROX index. *Journal of Critical Care*. 2016 Oct;35:200–5.

Dive deeper. Exclusive offers for HFNC enthusiasts

www.hamilton-medical.com/HFNC-upgrade

Learn more



Imprint

Published by: Hamilton Medical

Author: Kathrin Seeger

Release date: August 2022

Edition: 2



More information:

www.hamilton-medical.com/hiflow



Manufacturer:

Hamilton Medical AG

Via Crusch 8, 7402 Bonaduz, Switzerland

+41 (0)58 610 10 20

info@hamilton-medical.com

www.hamilton-medical.com

ELO20150701N.02

The products shown here are not available for purchase by the general public. The information provided here is intended for healthcare professionals only. Always read the labels and follow the product's instructions for use. Specifications are subject to change without notice. Some features are options. Not all features are available in all markets. All images are used for illustrative purposes only and may not accurately represent the product or its use. For all proprietary trademarks (®) and third-party trademarks used by Hamilton Medical AG see www.hamilton-medical.com/trademarks. ©2025 Hamilton Medical AG. All rights reserved.

The basics of high flow nasal cannula therapy